


HPLC Analysis of Dehydroascorbic Acid and Ascorbic Acid (Vitamin C) on SeQuant ZIC-HILIC

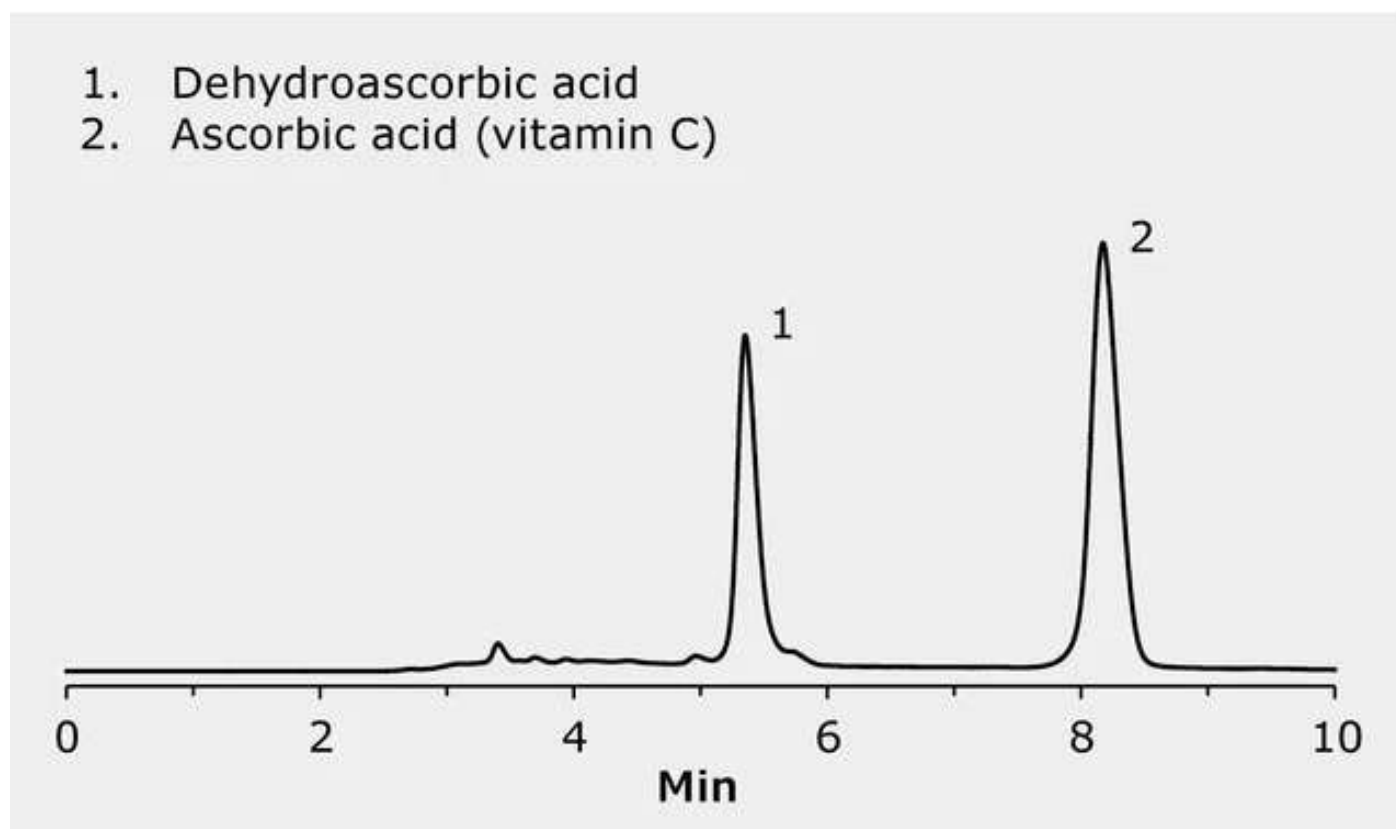

Conditions

column	SeQuant ZIC-HILIC, PEEK, 150 x 4.6 mm, 5 µm, 200 Å... (1.50455.0001)
mobile phase	[A] acetonitrile; [B] ammonium acetate; (70:30, A,B)
flow rate	0.5 mL/min
pressure	375 psi (2.6 MPa)

Conditions

injection	5 µL
detector	240 nm

Description

Featured Industry	Vitamins, Nutraceuticals, and Natural Products
suitability	application for HPLC